


# TENKER

MINI MASSAGER SM9126N

## USER MANUAL

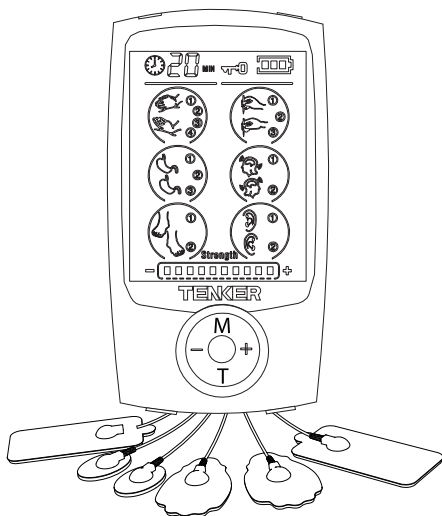

**Explore Your Adventure**

[www.tenker.co/warranty](http://www.tenker.co/warranty)

[facebook.com/Tenkerofficial](https://facebook.com/Tenkerofficial)

## ACTIVATE YOUR 12 MONTH WARRANTY & GET EXCLUSIVE GIFT

Register within 2 weeks after receiving your new product.

**PLEASE VISIT NOW**

→ [www.tenker.co/warranty](http://www.tenker.co/warranty)

**CUSTOMER SERVICE**

✉ [support@tenker.co](mailto:support@tenker.co)

📘 @Tenkerofficial



### **NOTE**

1. Avoid touching the pads when the device is turned on.
2. Do not place pads near the heart, head or carotid.
3. Not suitable for use by pregnant women, people with implanted devices or heart disease.
4. Apply a few drops of water when the pads have lost stickiness.

## A Big Thank You

We're thrilled that you have chosen to purchase the Electronic Massager Machine. At TENKER, we're committed to providing you with products that enhance your daily life. We also aim to give our customers excellent support. Should you have any questions or feedback about your new product, simply get in touch with our customer support team at 1-855-558-8666, or by email at [support@tenker.co](mailto:support@tenker.co).

Enjoy your new product!



# CONTENTS

Introduction .....	01
Safety Warning .....	02
Transport and Storage Symbols .....	03
Cleaning and Maintenance .....	04
Product Specifications .....	04
Set up .....	05
Operating Instruction .....	08
Mode Instruction .....	09
Recommended Use Positions.....	10
Recommended Practice .....	12
Trouble Shooting .....	13
Warranty Card .....	15
Contact Information .....	16

## INTRODUCTION

This Electronic Pulse Massager is a user friendly product. This Electronic Pulse Massager delivers electronic impulses to tired and sore muscles. These impulses are generated by the device and delivered through electrodes to the direct proximity of the target muscles. The impulses mimic the action potential coming from the central nervous system to trigger contraction of the muscles. Some programs will improve resistance to fatigue and increase endurance, while others can improve force produced by the muscle.

This device is portable, compact and stylish, as well as being controlled by a micro-processor. It may be helpful in relieving minor aches and pains in various parts of the body such as the waist, shoulders, Joints, hands and feet. The device is especially helpful in relieving stiff neck muscles and sore muscles from over exercising. Using it responsibly and properly will provide you with a relaxing pulse massage. You will feel a real difference this device can make to your life.

## **SAFETY WARNINGS**

### **Danger**

1. Do not use this unit if you are using an implanted electronic medical device such as pacemaker or other life-maintaining device.
2. Do not use this unit if you are wearing electrode medical devices for monitoring health functions.

### **Warning**
















1. If you have any of the following conditions, please consult your physician before purchasing or using this device:  
Acute conditions, malignant tumors, infectious diseases, pregnancy, implanted devices, heart disease, high fever, abnormal blood pressure, lack of skin sensation or numbness in the extremities, abnormal skin conditions or any condition requiring the active supervision of a physician.
2. Please remove the battery from the device when it is not in use.

### **Notice**

1. Do not use this device while driving.
2. Do not use this device while sleeping.
3. Do not use this device in areas with high humidity such as a bathroom.
4. Keep this device out of reach of children.
5. Stop using this device and consult your physician at once if you feel pain, discomfort, dizziness or nausea.
6. Do not attempt to move the electrode heads while the device is in operation.
7. Do not use the device near the heart, on the head, near the mouth, pudendum or blemished areas of skin.
8. Keep the device away from wet, hot place and out of direct-sunlight.
9. Store the device in a cool, dry place.
10. Do dispose of batteries according to the local laws and regulations of your region.

11. Medical electrical equipment needs special precautions regarding
12. Portable and mobile RF communication equipment can affect medical electrical equipment.

## TRANSPORT AND STORAGE SYMBOLS

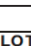	Fragile, handle with care		Type B applied part
	Keep the product in a dry place, away from water and rain.		Attention, see operating manual for user instructions.
	This way up		Manufacturer
	Product package should be recycled		Symbol for "AUTHORISED REPRESENTATIVE IN THE EUROPEAN COMMUNITY"
	Unrecyclable		CE marking, Certificate issued by SGS.
	Date of manufacture		FDA 510(K) approved
	Serial number		IP code of the device
	Batch code		



## **CLEANING AND MAINTENANCE**

### **Device**

1. Please use a lightly moistened cloth to clean the device and wipe gently.
2. Make sure the device is off when you clean it.
3. Do not let water to get into the device.

### **Pads**

1. Make sure both pads not overlapping and make sure your skin is free from any dirt, oil or lotion.
2. Applying a few drops water before and after usage will increase the longevity of the pads.
3. Place pads on pad holder or place plastic film on the sticky side of pad to keep the gel moisture.

## **Product Specifications**

### **Included in this standard package:**

1. TENS device x 1 pc
2. Large pad x 2pcs
3. Medium pad x 4pcs
4. Small pad x 2pcs
5. Electric wire x 3pcs  
(when use 4 line cable, make sure to use at least 1 pc white(positive)  
+1 pc blue line (negative) at the time.)
6. USB cable & AC adapter x 1 set
7. Instruction manual x 1 pc
8. Collection board for pad and wire x 1 pc
9. Bag x 1 pc


## Features

- Big size blue LCD screen
- 2 outputs with 16 modes
- 20 level intensity
- 20 minutes default timer, can be adjustable from 10-60 minutes
- Rechargeable lithium battery
- Battery indicator to know when need to recharge battery
- AC adapter 110-240v, USB Charger 5v


## SET UP

Unpack the box of the product, take the product and accessories out, and connect the electrode pad onto the device.

1. The device has two independent channels A&B . Make sure which channels (channel A or B or both) you want to use before you connect wires to the device.
2. Put pads on your skin before turn on the device. Use at least one pair of pads, one pad alone will not work.

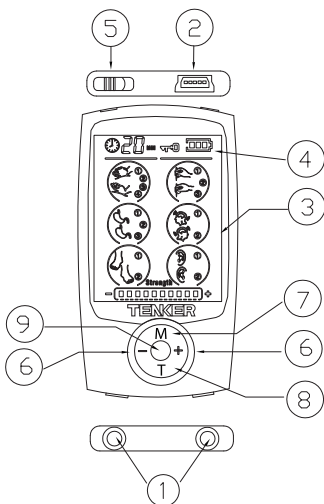

- |                           |                         |
|---------------------------|-------------------------|
| 1. Two output connector   | 2. USB port             |
| 3. LCD screen             | 4. Battery indicator    |
| 5. On/ off switch         | 6. Intensity adjustment |
| 7. Mode selector          | 8. Time adjustment      |
| 9. Built-in Mode Selector |                         |

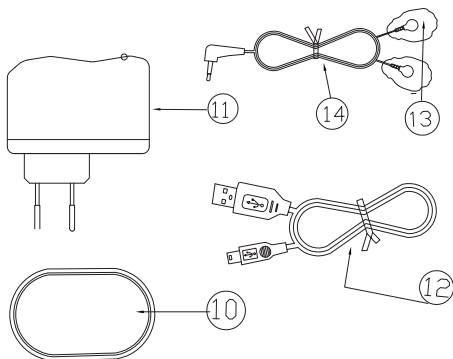

10. Pads- sticking board

12. USB cable

14. Pad wires

11. Charger

13. Massage pads

## OPERATING THE PRODUCT

1. Connect the pads to the wire by clipping them on, and then remove the protective film.
2. Place the pads where the massage sensation is desired, and plug into one of the ports on the device. For additional relief, repeat the procedure by putting another two pads on the desired area and plug into the second port.

**Make sure both pads on the skin and not overlapping and make sure your skin is free from any dirt, oil or lotion.**

3. Turn the device on by sliding the on/off switch at the top of the unit to the 'on' position. The unit will be now activated.
4. Press the "M" button to select a different main mode (there are 6 main modes to choose from, each main mode contains different sub-modes). Press the center button to select a different sub-mode within the current main mode.

**Please remember that the strength setting will default to the lowest setting with each mode selection you make.**

5. Increase the intensity by pressing the + button; decrease the intensity by pressing the - button.
6. Press T button, 10,20,30,40,50,60mins timers could be selected.

**Default time is 20 minutes.**

7. Turns off the device.
8. Remove the pads from your skin.

**Gently remove the pads from your skin and place them back on the provided pads holder or protective plastic film.**

**Gently remove by pulling the corner of the pads. Do not pull the leads.**

## MODE INTRODUCTION

There are several modes to choose from; The "M" cycles through the 6 main modes, while the silver middle button allows choices (from 2 to 4 varieties) within that mode.

There are 16 small modes in total:

### 1.Kneading

- 1)Lateral palm kneading
- 2)Flat palm kneading
- 3)Crosscut fist kneading
- 4)Vertical fist kneading

### 2.Acupuncture

- 5)Acupuncture
- 6)Moxibustion
- 7)Acupuncture kneading

### 3.Cupping

- 8)Rubbing reduction
- 9)Cupping
- 10)Fibrillation massage

### 4.Tapping

- 11)Lifting-pinching manipulation
- 12)Vibrating manipulation


### 5.Shiatsu-Deep

- 13)Thai massage
- 14)Shiatsu

### 6.Auto-Combination

- 15)Auto-Combination Massage Mode 1
- 16)Auto-Combination Massage Mode 2

## RECOMMEND USE POSITIONS

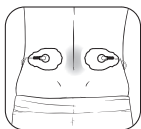

### LOWER BACK

Mode 1,4,5,6,7,10,11

Attach both pads on the lower back with the backbone in the center.  
Do not place on the backbone or spine.



### SCIATIC

Mode 1,2,3,4,5

Attach on pad below and above the region in pain, both on the same side (NOT on spine)



### DELTOID

Mode 2,3,4


Attach one pad in front and one in back of the muscle.



### ARM


Mode 1,2,3,6,7,15,16

Attach both pads on either side of the area.

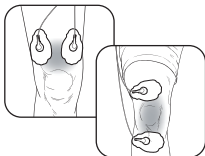


### KNEE

Mode 1,2,3,4,5,6,7,10  
Attach both pads above the knee  
or above and below joint.

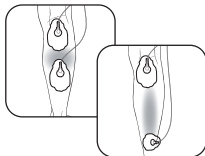


### CALF

Mode 1,2,3,4,5,6,7,8,14,15  
Attach both pads on the calf/leg  
where you feel pain.  
(Do not place electrode pads simultaneously  
to the calves of both legs.)

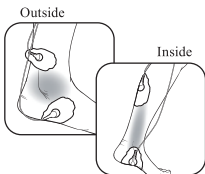

### ANKLE/FOOT

Mode 1,2,14,15  
Attach pads per the illustration, on the left/  
on the outside of your ankle/foot.  
Attach the pads per illustration on the right/  
on the other side of your ankle/foot

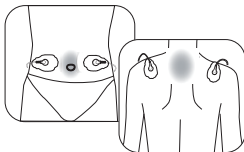

### ABS/SHOULDER

Mode 1,4,5,6,7,10,11,12,14,15,16  
Attach pads per the illustration.  
(Do not put pads on the spine)



## RECOMMENDED PRACTICE

- Duration suggested for each skin area is 20 min and 2 times per day. Consult with your physician for longer and more frequent uses.
- Start from the lowest intensity and gradually adjust the intensity to a comfortable level at a scale from 1 to 20.
- Good skin care is important for a comfortable use of device. Be sure the treatment site is clean of dirt and body lotion.
- Keeping the electrode in the collector after use will extend its lifespan. The electrode is disposable and should be replaced when it loses the adhesiveness. To purchase additional electrodes, please contact the seller.

## TROUBLE SHOOTING

If your device is not operating properly, please check below for common problems and suggested solutions.

If the recommended action does not solve the problem, please contact the seller at

Official email : support@tenker.co

Tel No: 1-855-558-8666

Problem	Possible Cause	Solution
One pad feels stronger than the other	This is normal. Different area of your body will react differently	Nothing needs to be done. Make sure the pads are moist and making good contact.
The intensity is not felt with a very weak intensity level	Pads are not attached to the body firmly	Attach both pads firmly to the skin
	The transparent films are still stuck to the pads	Peel off film on the adhesive surface of pads
	The pads stack together or overlap	Do not stack pads together or overlap pads
	The cord is not properly connected to the unit	Connect cord correctly into the jack
	The intensity setting is getting weak	Increase the intensity level
	The battery capacity is low	Charge the battery
The skin turns red or the skin feels irritated	The adhesive surface of the pads is dirty or dry	Wash adhesive surface of pads gently with your fingertips for about 3 seconds under slow running water
	The therapy time is too long or the intensity is set too high	Reduce the application time or reduce the intensity
	The electrode pad surface is worn out	Replace electrode pad

Problem	Possible Cause	Solution
No power source; no display on LCD.	The battery capacity is depleted	Charge the battery
Power cuts off during use	The battery is weak	Charge the battery
	The cord is broken	Replace the cord
It is difficult to attach the pad to the skin	Have you remove the transparent film from the pad?	Peel off film on the adhesive surface of pads
	Was the pad applied immediately after washing?	Dry the pad
	Is the adhesive surface of the pad damaged?	Replace the pad
Adhesive surface of pad is not sticky	Pads get deteriorative	Contact the vendor for replacements.
	Were the pads stored under high temperature, high humidity, or direct sunshine?	Replace the pad.

## WARRANTY CARD

Certificate of Quality Warranty Card	
Product model	SM9126N
Product name	Electronic Mini Massager
Date of purchase	
Purchase Location	
Retailer's Name	
Retailer's TEL	
Problems of the Product	

## CONTACT INFORMATION



### Manufactured by:

Hong Qiangxing (Shenzhen) Electronics Limited Address:4F,  
Jingcheng Building,Xicheng Industrial Area, Xixiang  
Street,Bao'an District,Shenzhen,  
Guangdong,China

EC

REP

### Wellkang, Ltd

Suite B, 29 Harley Street London, United Kingdom W1G 9QR  
OMGD INC  
ADD:245 East Main Street,Suite 115,Alhambra,CA 91801

Made in china

**WEEE-Nr : DE83395812**

## **ACTIVATE YOUR 12 MONTH WARRANTY & GET EXCLUSIVE GIFT**

Register within 2 weeks after receiving your new product.

**PLEASE VISIT NOW**

→ [www.tenker.co/warranty](http://www.tenker.co/warranty)

**CUSTOMER SERVICE**

✉ [support@tenker.co](mailto:support@tenker.co)

📘 @Tenkerofficial



### **NOTE**

1. Avoid touching the pads when the device is turned on.
2. Do not place pads near the heart, head or carotid.
3. Not suitable for use by pregnant women, people with implanted devices or heart disease.
4. Apply a few drops of water when the pads have lost stickiness.



# TENKER

✉ support@tenker.co

🖱 www.tenker.co

☎ 1-855-558-8666 (Mon - Fri 9:00 AM - 5:00 PM EST)

WEEE-Nr: DE83395812

CE 0123



Made in China